



Woodruff Scientific, Inc

4000 Aurora Ave N, Suites 5 & 6,

Seattle, WA 98103

(206) 905 9477 8am to 5pm Pacific


sales@woodruffscientific.com

<http://www.woodruffscientific.com>

Model number(s): R1-F

Descriptive name: Fiber-Coupled Interferometer

Optical setup of a 2-chord, fiber-coupled interferometer on a 12x24 in optical table (enclosure not shown)

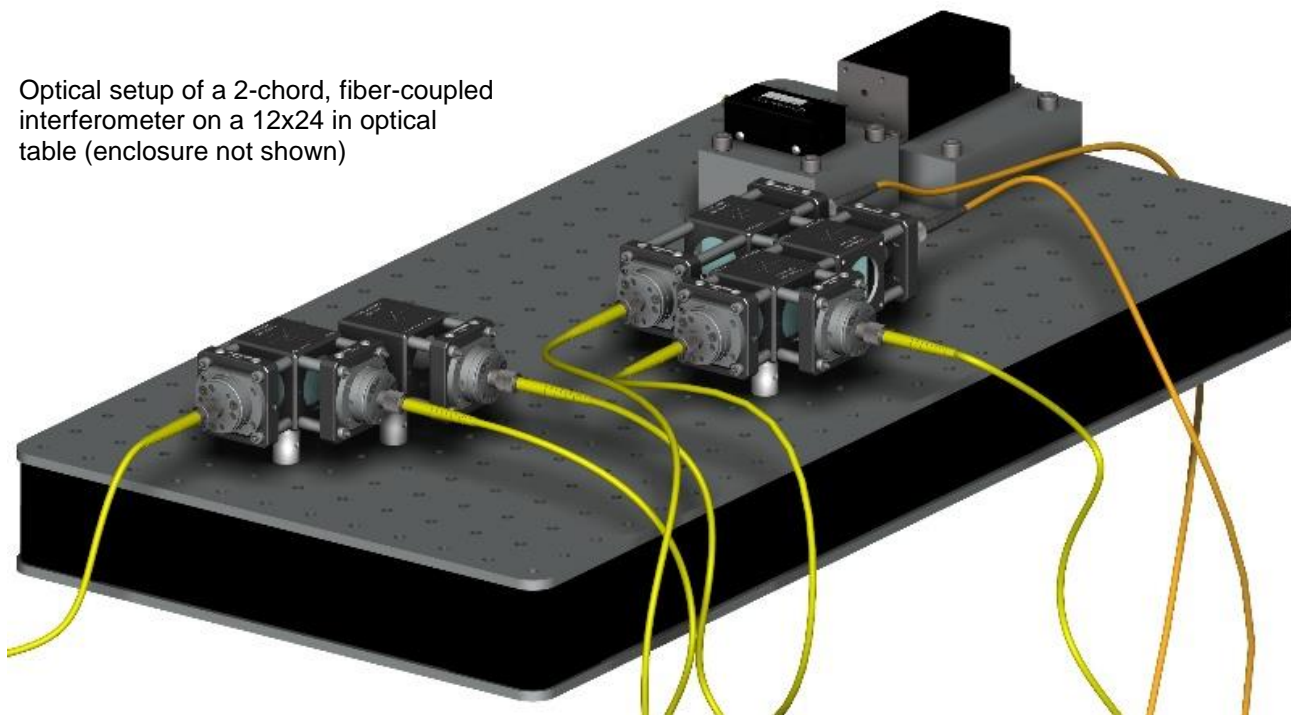

Features:

- Measures chord-averaged electron density
- Fiber-based design allows line of sight to be changed without realignment of optics
- High-coherence laser allows for single reference beam and large path length discrepancies
- High-frequency laser provides high temporal resolution
- Heterodyne configuration for extended unambiguous measurement range
- Extendable to multiple chords
- Laser sources in 532 nm – 5 μ m wavelengths
- Uses FC/PC fiber optic connections for high quality beam launching
- Safe operation with interlocked/key-switched enclosure and low-power setting for beam alignment

© Woodruff Scientific Inc, 4000 Aurora Ave N, Suites 5&6, Seattle, WA 98103

Sales: sales@woodruffscientific.com (206) 905 9477



Woodruff Scientific, Inc

4000 Aurora Ave N, Suites 5 & 6,
Seattle, WA 98103
(206) 905 9477 8am to 5pm Pacific

sales@woodruffscientific.com

<http://www.woodruffscientific.com>

Model number(s): R1-F

Descriptive name: Fiber-Coupled Interferometer

Operational ratings:

| | 0.532 μm | 5 μm |
|--|---|--|
| Electron density (n_e): | $< 4 \times 10^{27} \text{ m}^{-3}$ | $< 5 \times 10^{25} \text{ m}^{-3}$ |
| Density resolution* ($\delta\langle n_e L \rangle$): | $1 \times 10^{19} \text{ m}^{-2}$ | $1 \times 10^{18} \text{ m}^{-2}$ |
| Temporal resolution** ($d\langle n_e L \rangle/dt$): | $3 \times 10^{30} \text{ m}^{-2}\text{s}^{-1}$ | $3 \times 10^{29} \text{ m}^{-2}\text{s}^{-1}$ |
| Path length mismatch: | $< 2 \text{ m}$ between reference and probe beams | |

* corresponding to 1.0° phase resolution; actual resolution may vary with application

** depending on data acquisition

Options:

- *Laser wavelength in the range 532 nm – 5 μm*
Higher frequencies allow greater density before reaching cutoff, faster temporal response, and better spatial resolution. Lower frequencies are have better density resolution, are less sensitive to mechanical motion, and can be more economical in many applications.
- *Number of chords*
Having multiple chords allows for the spatial reconstruction of the electron density profile or analysis of plasma propagation, but increases the cost and complexity of the system.

WSI also offers a 2-color, fiber-coupled interferometer for vibration compensation, as well as standard HeNe, CO₂, and microwave interferometers. [See the website](#) for more information.

[Contact WSI](#) for pricing, operational ratings specific to your application, and more information on which options may be appropriate for you. Please provide any information you may have about your application's density, size, and timescale, as well as any preferences.




Woodruff Scientific, Inc

4000 Aurora Ave N, Suites 5 & 6,
Seattle, WA 98103
(206) 905 9477 8am to 5pm Pacific
sales@woodruffscientific.com
<http://www.woodruffscientific.com>

Model number(s): R1-F

Descriptive name: Fiber-Coupled Interferometer

Schematic (2-chord option shown):

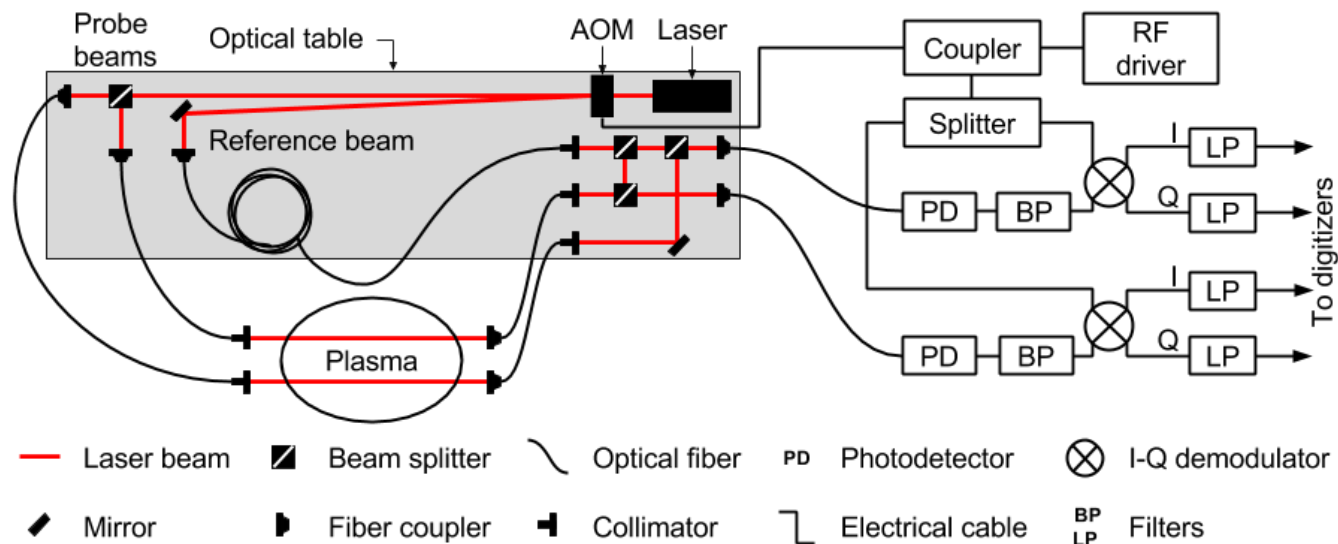

Shown above are the optical (left) and electrical (right) schematics for a 2-chord R1-F interferometer. The laser is sent through a standard AOM and split into the reference beam (now at different frequency) and the probe beam. The probe beam is then split into the desired number of chords. All beams are attenuated to the same power and coupled into fiber-optic cables. The cables are run to the desired viewing locations on the chamber, where a collimator is used to launch the probe beams into the plasma. The probe beams are collected on the other side, coupled back into fiber optic cables, and returned to the optics assembly. The reference beam is split and recombined with the probe beams. The recombined beams are converted to an electrical signal, filtered, and mixed with a local oscillator signal. The I and Q for each probe beam are output to coaxial cables for connection to the experiment DAQ system.



Woodruff Scientific, Inc

4000 Aurora Ave N, Suites 5 & 6,

Seattle, WA 98103

(206) 905 9477 8am to 5pm Pacific


sales@woodruffscientific.com

<http://www.woodruffscientific.com>

Model number(s): R1-F

Descriptive name: Fiber-Coupled Interferometer

Example data:

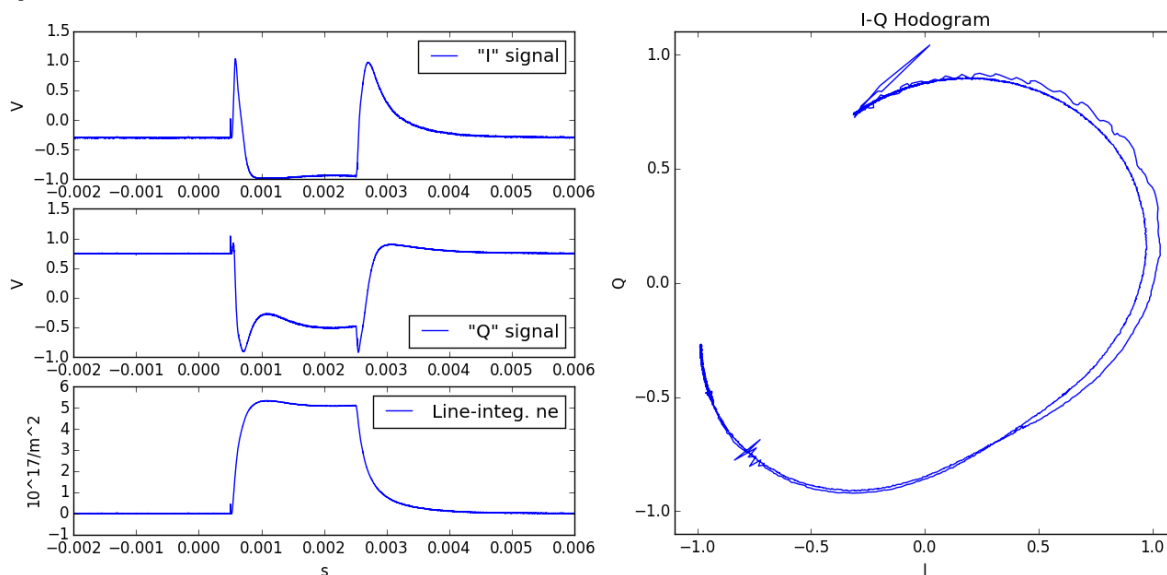

At left is example data from an interferometer. The I and Q signals represent the cosine and sine of the phase shift, respectively. The plot of Q vs. I demonstrates this relationship; ideal signals would trace a constant-radius arc or circle, the real-world plasma distorts the radius of the arc on the hodogram without changing the polar angle. The line-integrated electron density is shown bottom left.